


Remote alarm units for medical gases compatible with HGM devices



On-plaster alarm unit

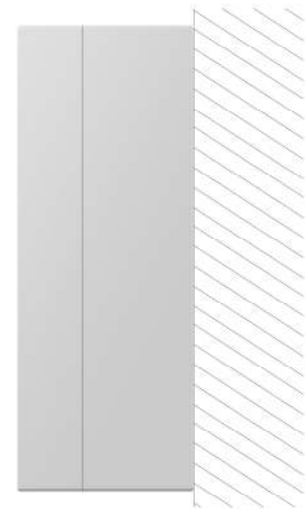

Dimensions (WxHxD)
171x121x55 (mm)



Under-plaster alarm unit



Dimensions (WxHxD)
174x137x5 (mm)

Function:

Reads from the master alarm unit all shared informations about medical gases. It shows exact copy of screen of master unit which is connected to. It has visual and audio alarm as same as alarm unit mounted in valve box. ALL devices are connected with use of FTP cables and RJ-45 plugs.

SYSTEM ADVANTAGES:

- limitless possibilities of expansion
- safe electrical supply - 12V DC
- power supply can be delivered to only one device and the whole network will be supplied
- one commonly available plug system RJ-45
- means to connect whole monitoring system to multiple central alarm units and/or PC and/or laptop
- compatible with 3rd party valve boxes with pressure switches or pressure sensors